

# Deconvolution in Scanning Electrochemical Microscopy

Dr. András Kiss

Department of General and Physical Chemistry  
University of Pécs, Hungary



Center for Integrative Physiology  
and Molecular Medicine  
Universität des Saarlandes, Homburg, Germany



Center for Integrative Physiology  
and Molecular Medicine


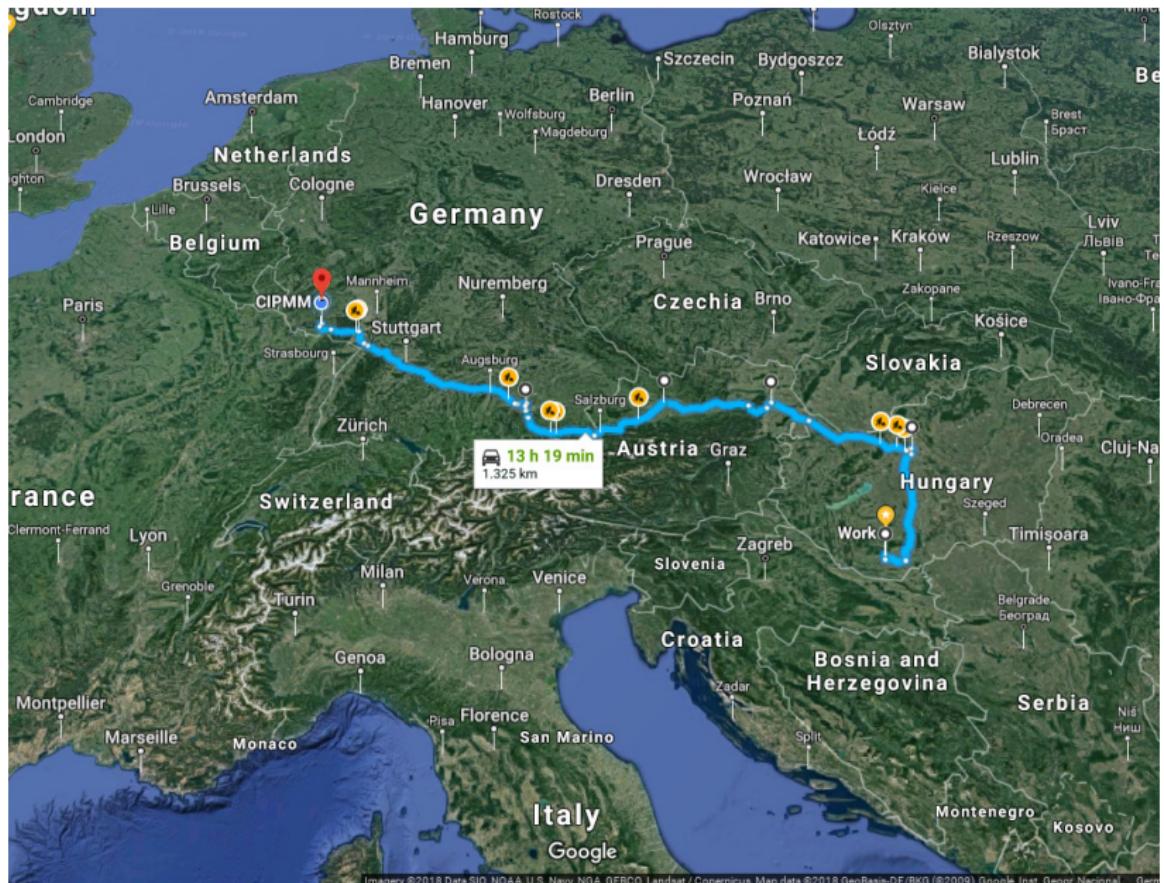



Image courtesy of NOAA/NESDIS/NGDC. Last updated: 2018-06-20T12:00:00Z

# *University of Pécs, Hungary General and Physical Chemistry*

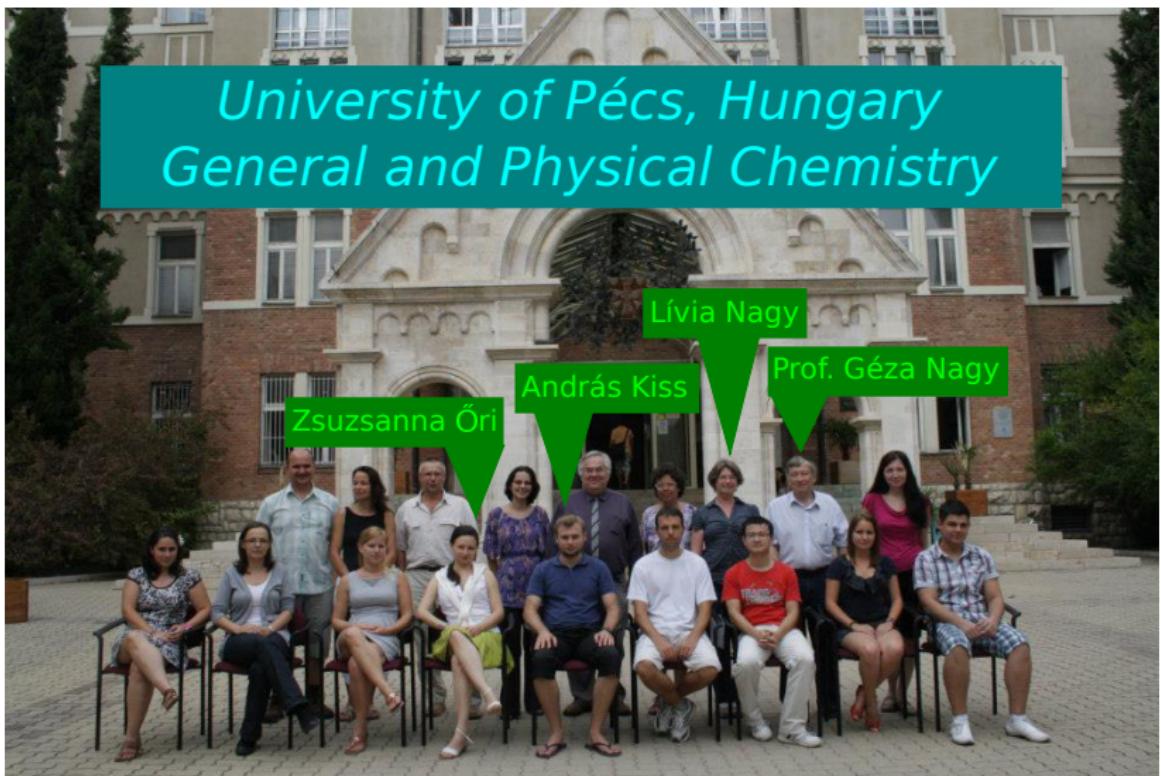


# Potentiometric Scanning Electrochemical Microscopy

A Scanning Probe Microscopic technique


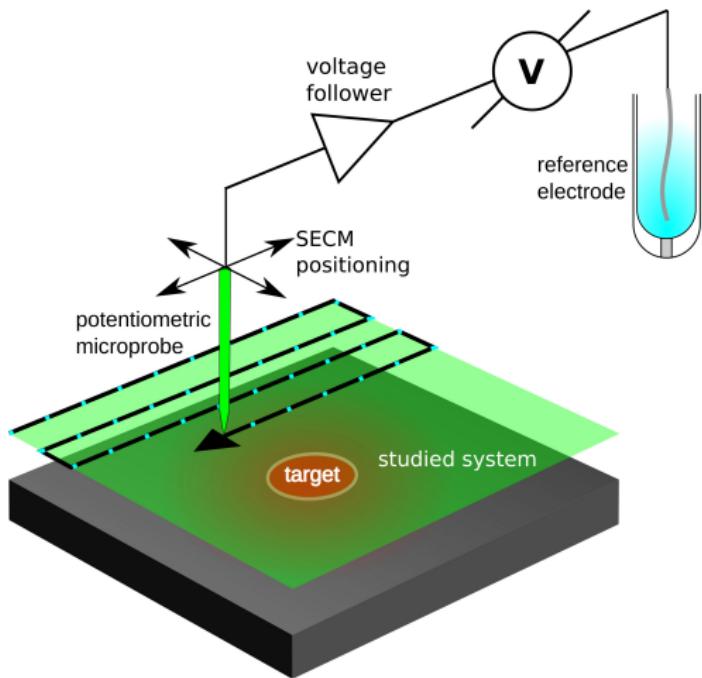


# Difference between conventional and scanning microscopic techniques

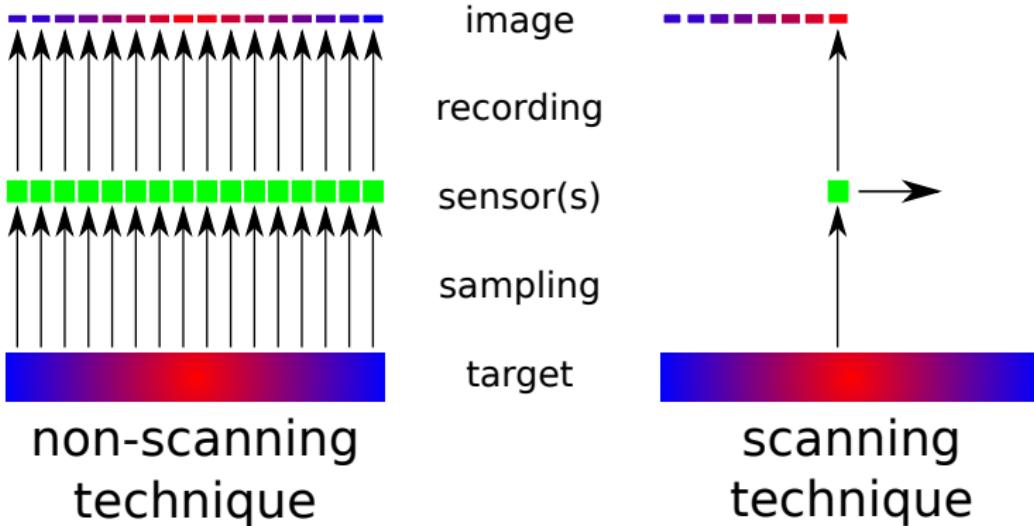

# Ion-selective micropipettes

As SECM probes

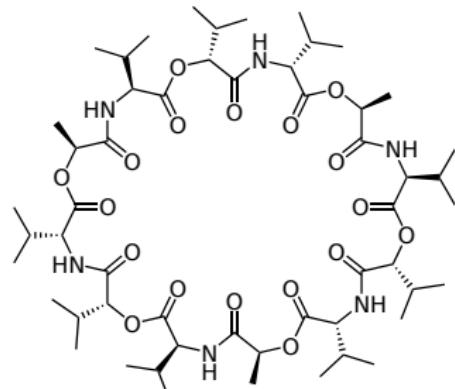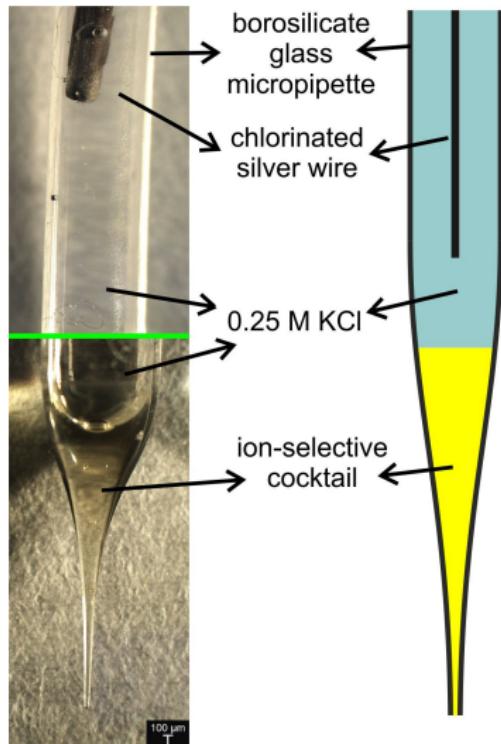

Valinomycin


$$E = E^\theta + \frac{RT}{z_i F} \ln \left[ a_i + \sum_j \left( k_{ij} a_j^{z_i/z_j} \right) \right]$$

Nikolsky-equation


# The problem with potentiometric SECM

Distortion at high scan rate

Slow

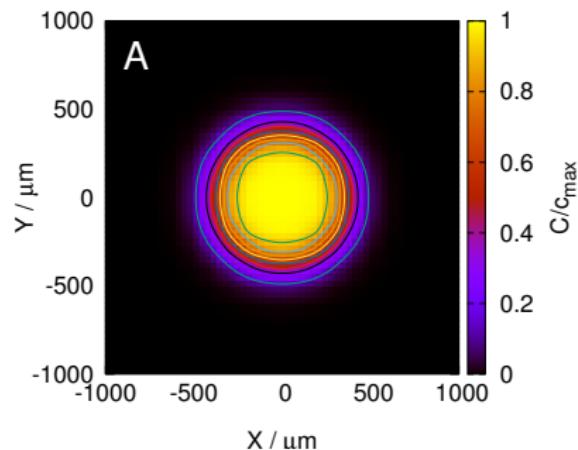


Fast

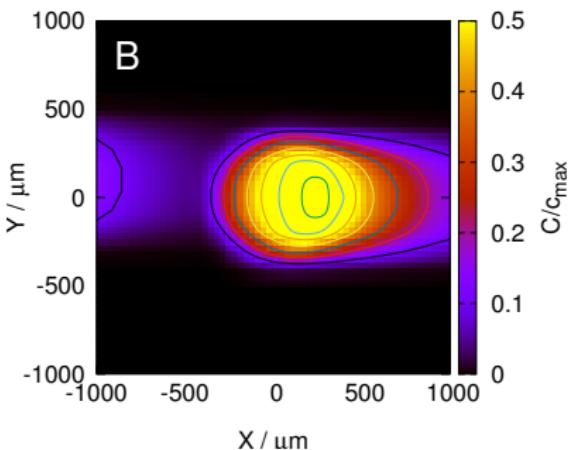


# Why is the image distorted?

Possible contributors to the lag

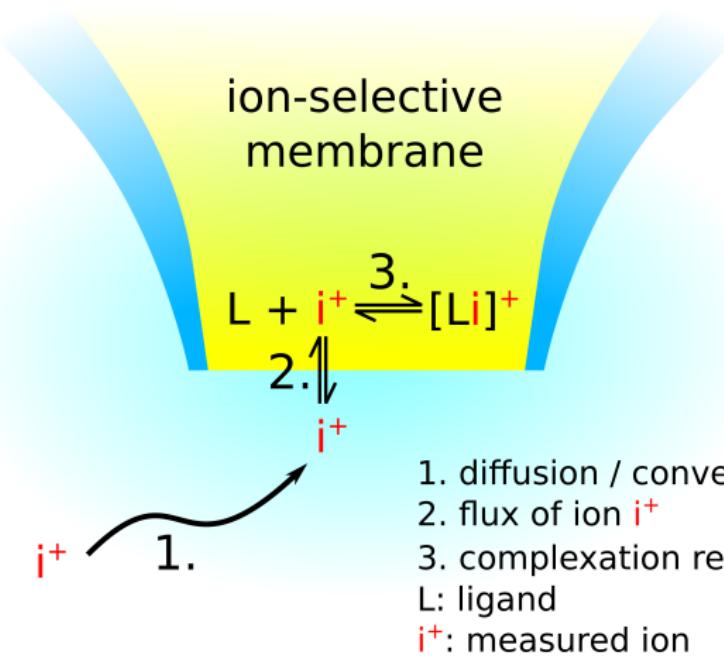

# Why is the image distorted?


The RC time constant

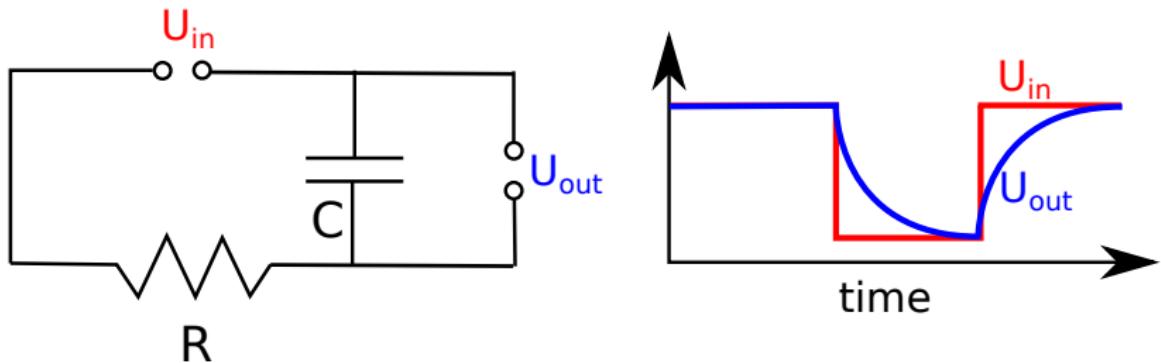

The time that is required to charge the capacitor by  $\approx 63\% (1 - 1/e)$ .


# Distortion of potentiometric imaging

In the case of a linescan

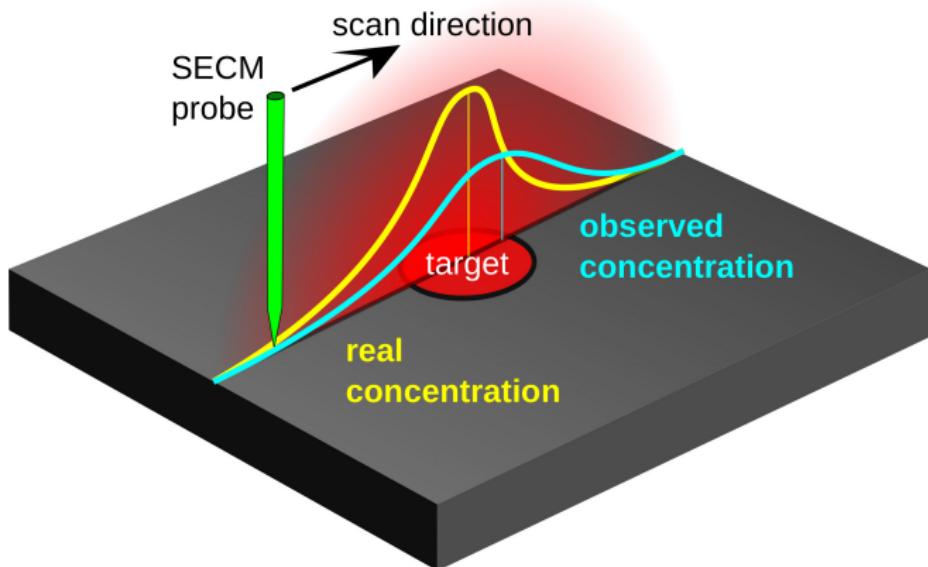


# Trade-off triangle of potentiometric SECM


Compromise between the three desired competing properties

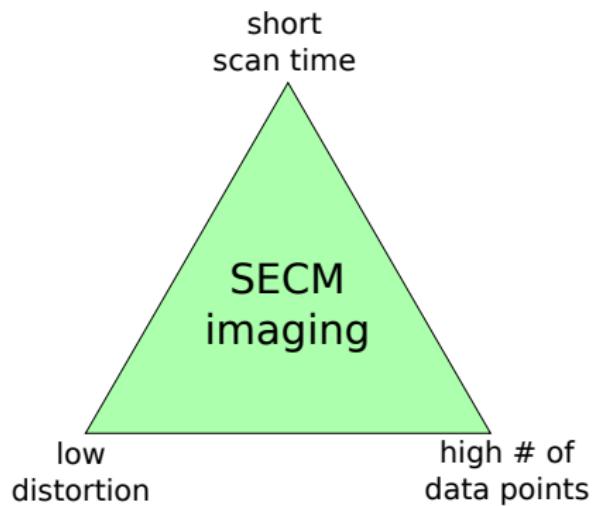

## The convolution function of the distortion

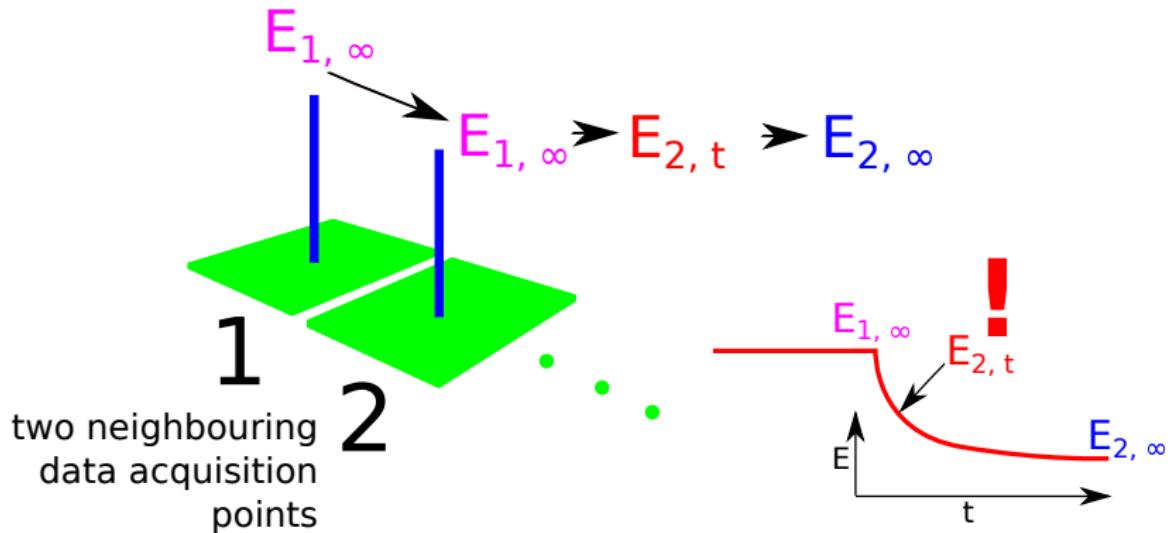


# Convolution and deconvolution

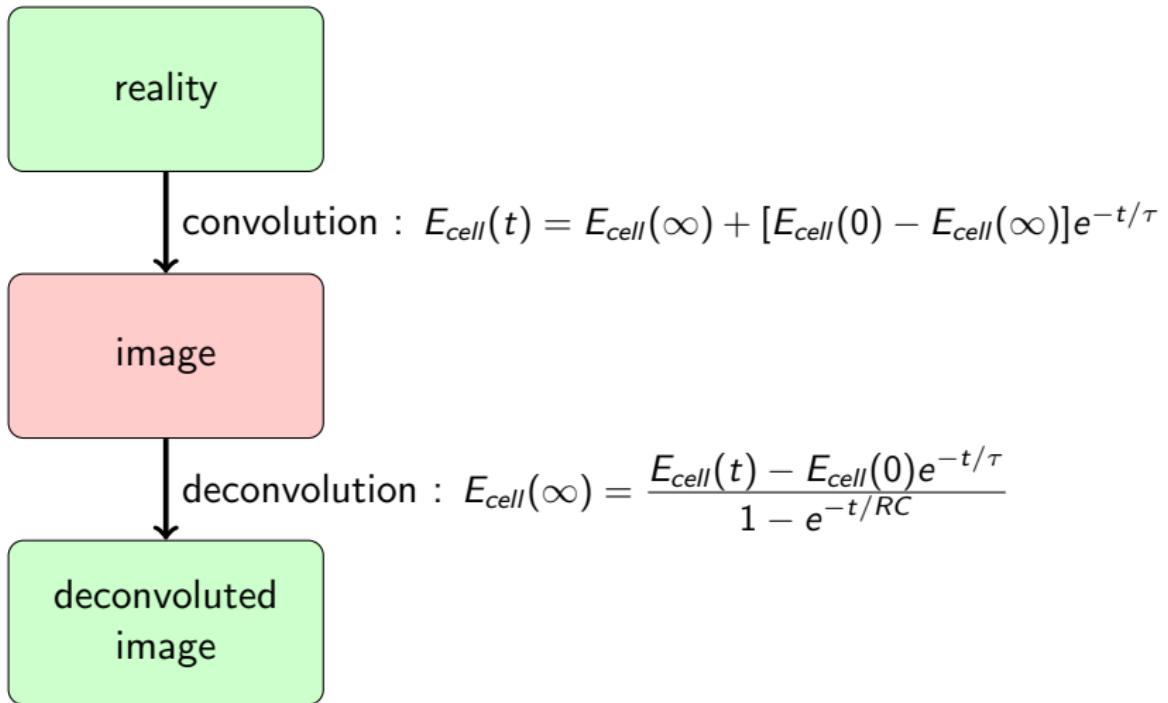


# Deconvolution of potentiometric SECM images

Recorded using the antimony microelectrode following the meander algorithm

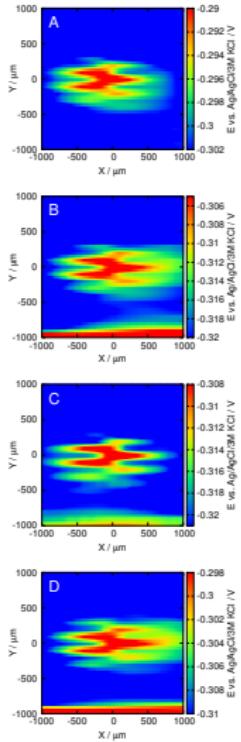


deconvolution  
→

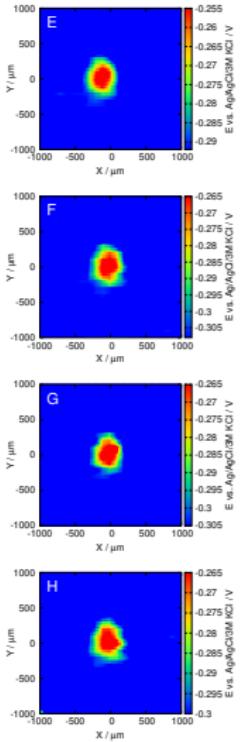

# Practical example: corroding carbon steel sample

Scanned with an antimony microelectrode

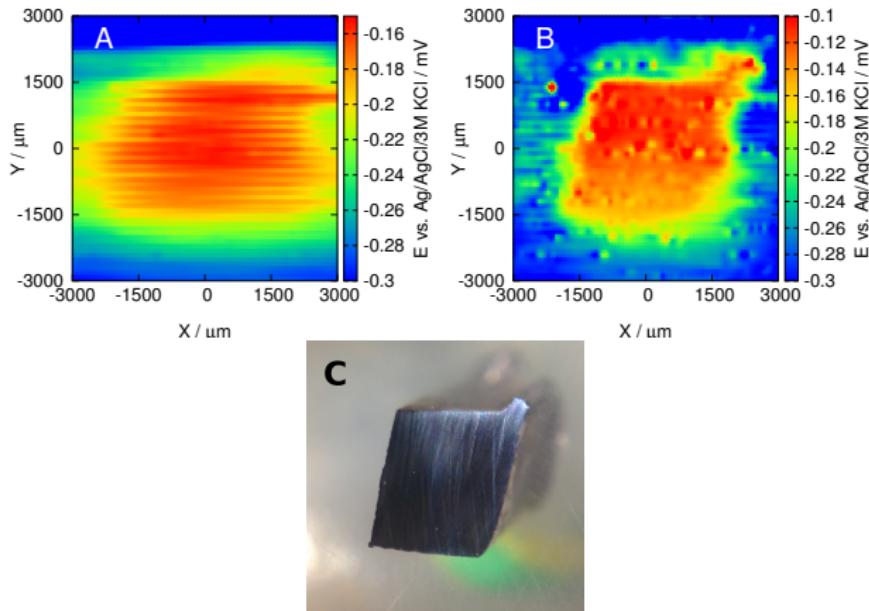

# 9th Workshop on Scanning Electrochemical Microscopy and Related Techniques

Warsaw, Poland, August 13-17, 2017.




Can it be done with amperometric SECM images?

# DAAD

Deutscher Akademischer Austausch Dienst  
German Academic Exchange Service

# Why would you want to do SECM imaging of a cell?

It's an easy way to make selective, non-invasive, high resolution, single cell experiments

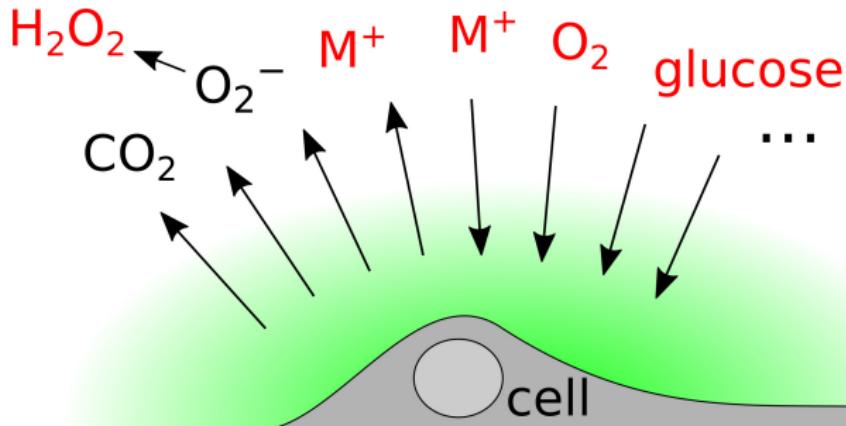

ANTIOXIDANTS & REDOX SIGNALING  
Volume 29, Number 6, 2018  
Mary Ann Liebert, Inc.  
DOI: 10.1089/ars.2016.6840


FORUM ORIGINAL RESEARCH COMMUNICATION

## Electrochemical Quantification of Extracellular Local $\text{H}_2\text{O}_2$ Kinetics Originating from Single Cells


Monika Bozem,<sup>1</sup> Philip Knapp,<sup>1</sup> Valentin Mirăceski,<sup>2</sup> Ewa J. Słowiak,<sup>1</sup> Ivan Bogeski,<sup>1,3</sup> Reinhard Kappi,<sup>1</sup> Christian Heinemann,<sup>4</sup> and Markus Hoth<sup>1</sup>





# The amperometric measuring cell

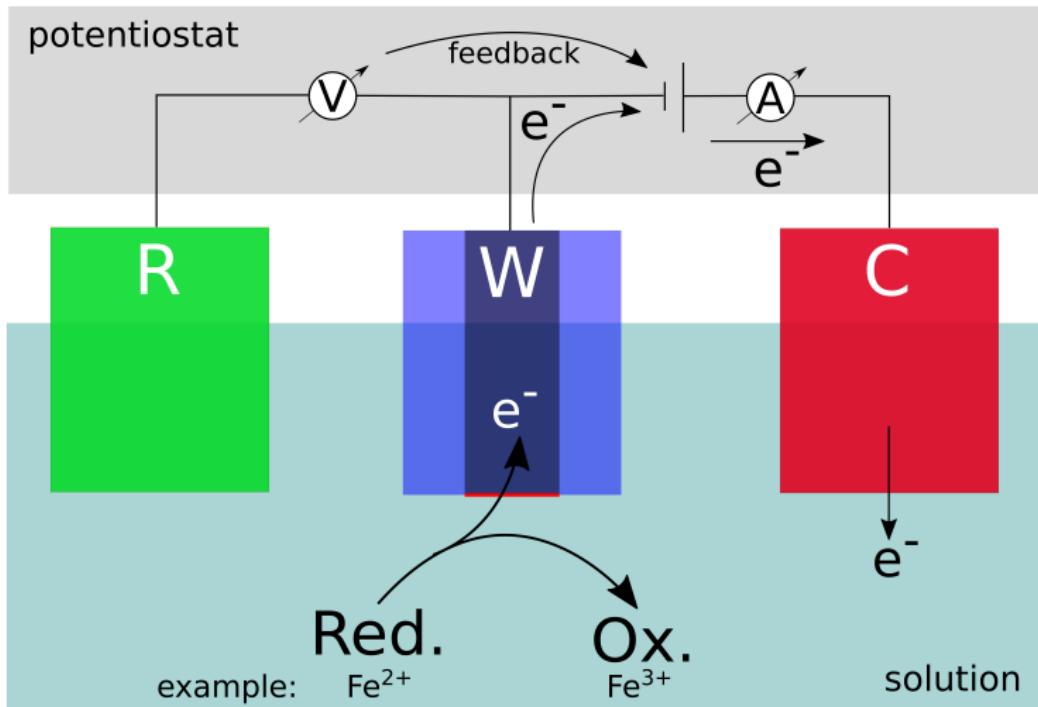

# Amperometric transient response

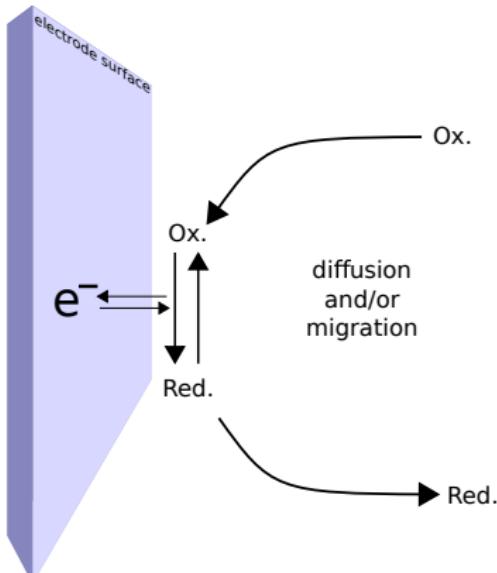


Simplified equivalent amperometric circuit

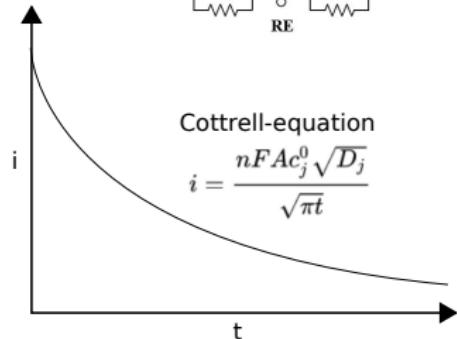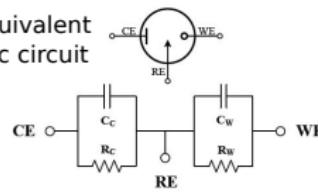


# Deconvolution of an amperometric image.

$\text{H}_2\text{O}_2$  oxidation current above a monocyte exposed to extracellular  $\text{H}_2\text{O}_2$  (10  $\mu\text{M}$ ).

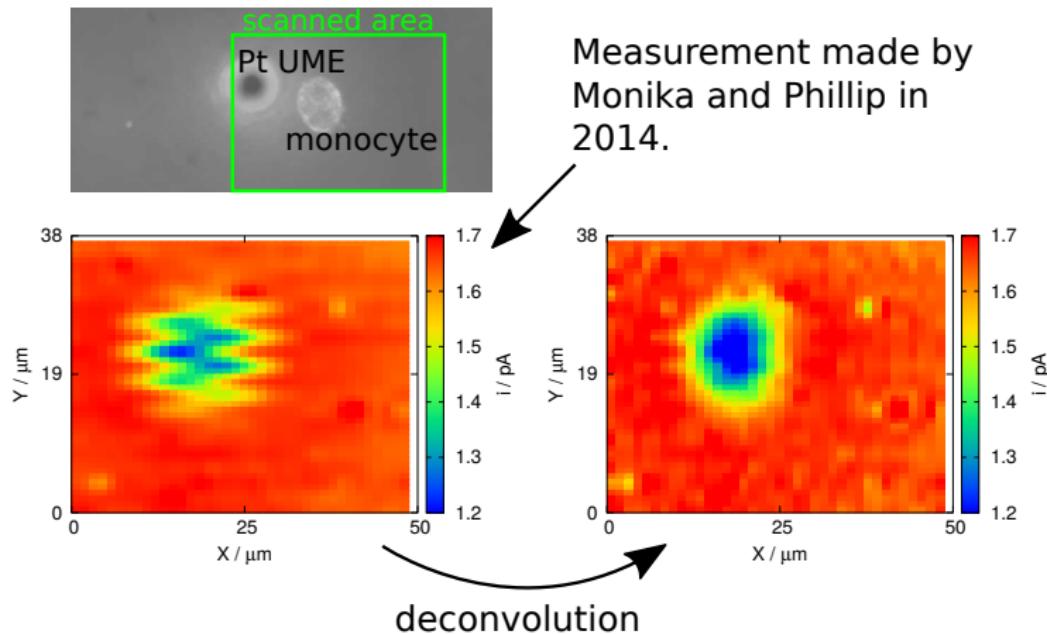


# Investigating a step response over a glass sheet edge.

With the ferrocene/ferrocenium system.

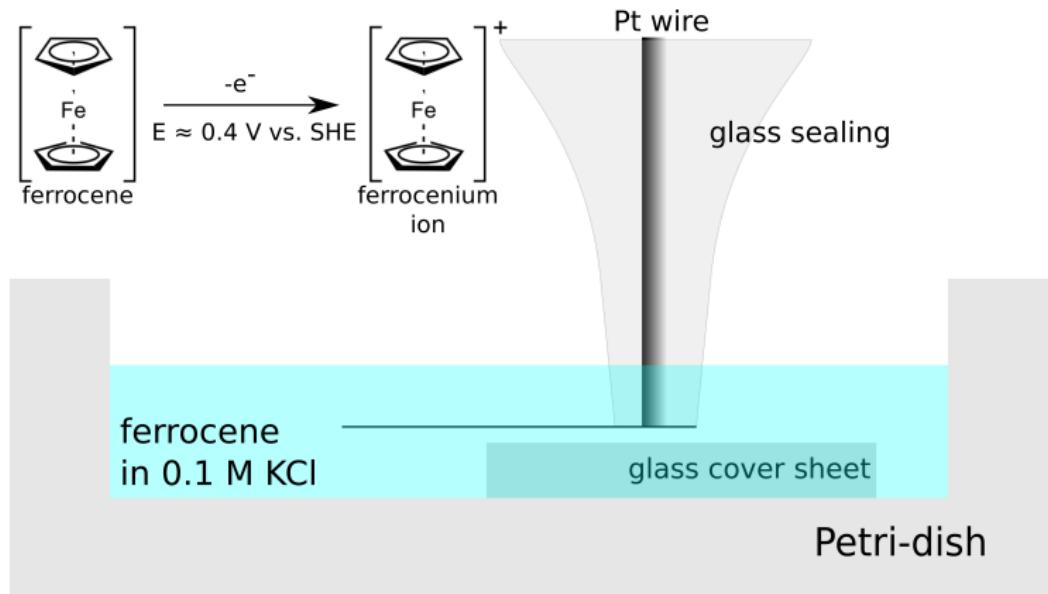


# Investigating a step response over a glass sheet edge.

## Results.

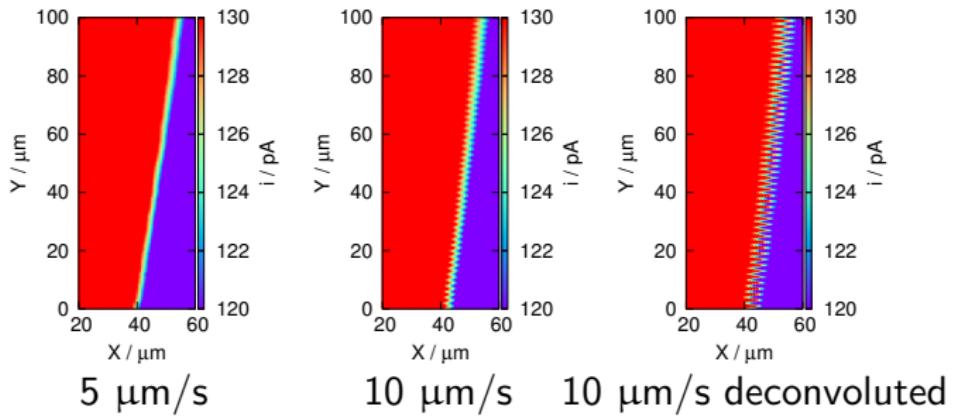


# Investigating a step response over a glass sheet edge.

Discrepancy caused by convection.

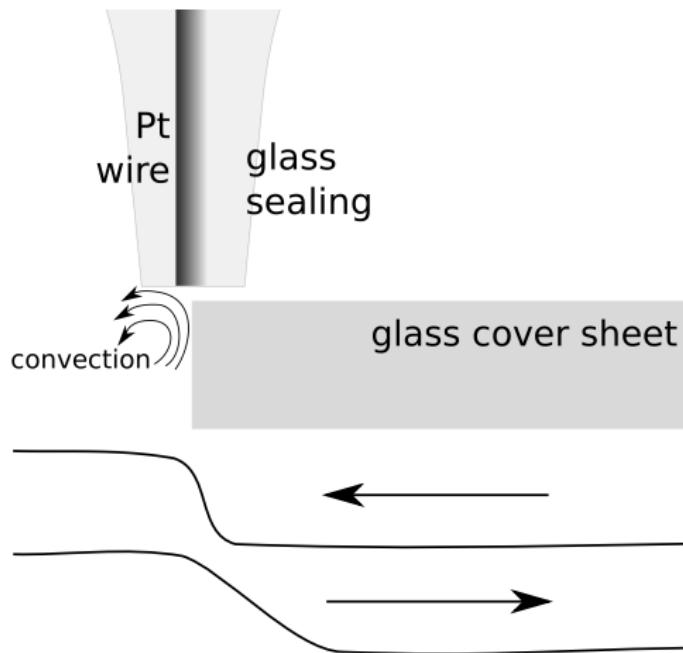



# Investigating a step response over a $d = 10 \mu\text{m}$ Pt wire.

With the ferrocene/ferrocenium system.

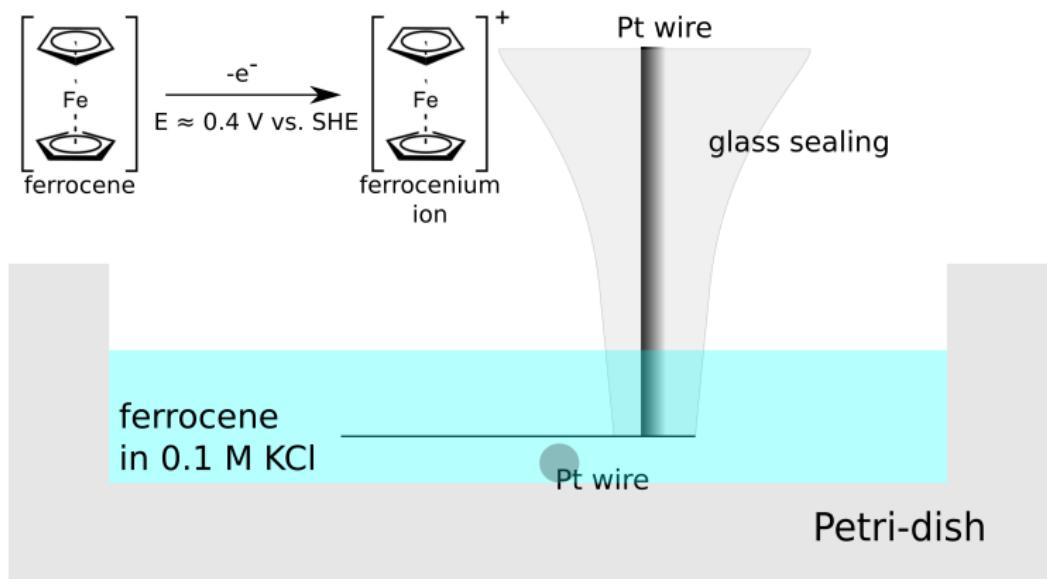


# Investigating a step response over a $d = 10 \mu\text{m}$ Pt wire.

Results.

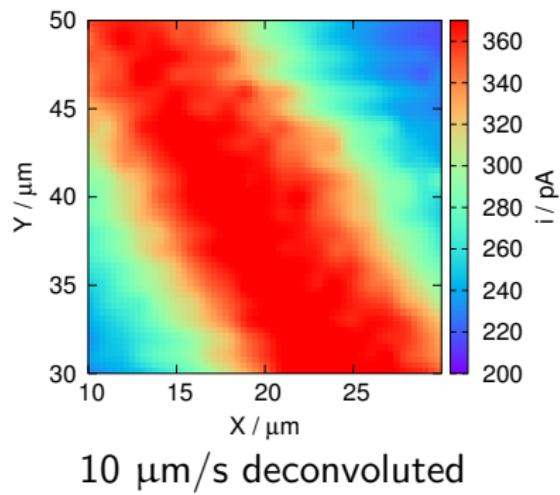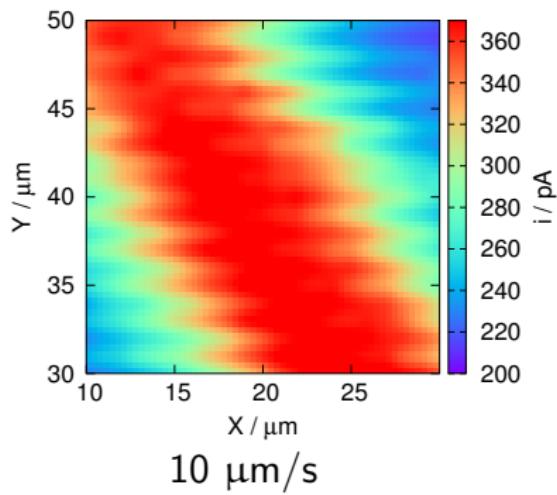

# $\text{H}_2\text{O}_2$ measurement over a TPA stimulated monocyte.



## Results.

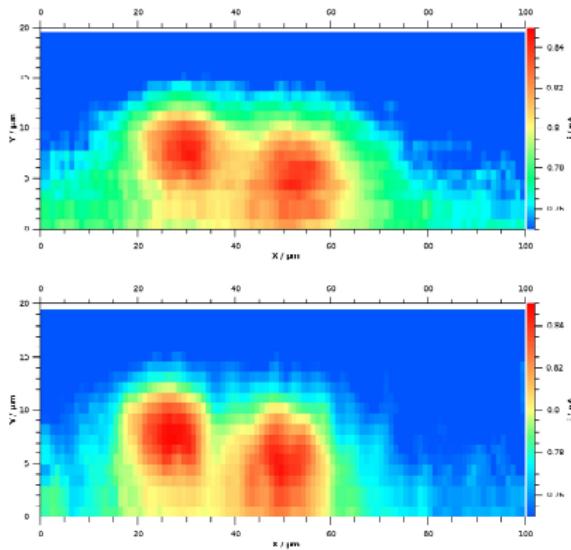

Top image: original measurement, made by Monika in 2014.  
Scanning speed was 2  $\mu\text{m}/\text{s}$ . Bottom image: deconvoluted image.



# Oxygen reduction current above a monocyte.

Experimental setup.

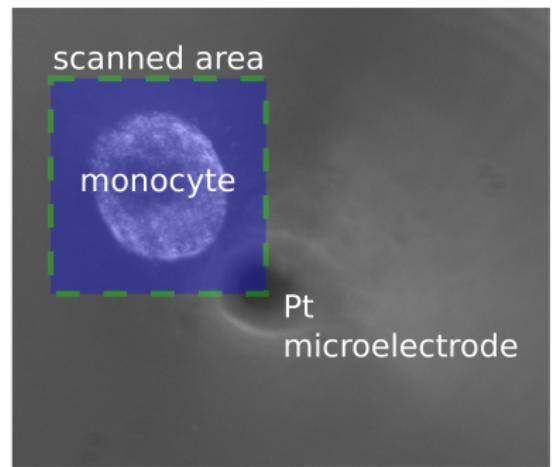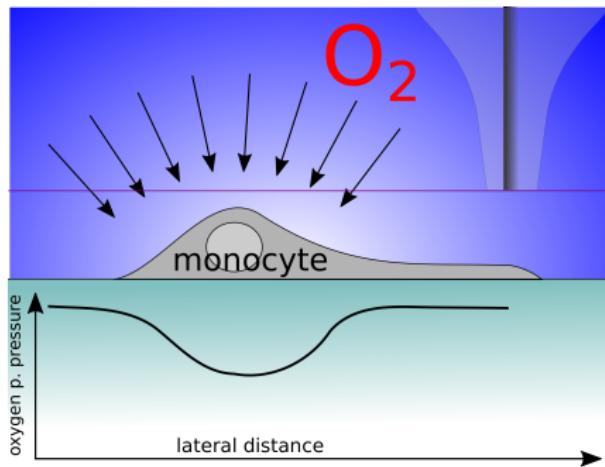

# Oxygen reduction current above a monocyte.

Results.

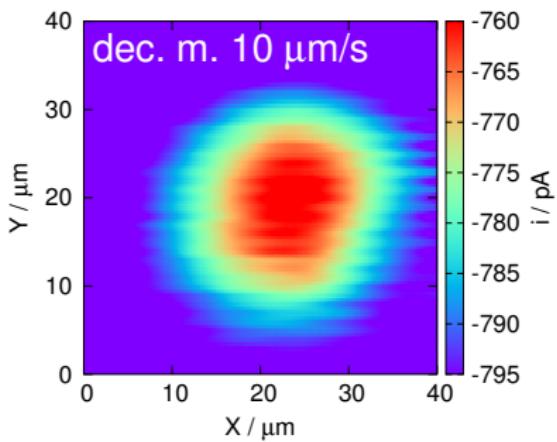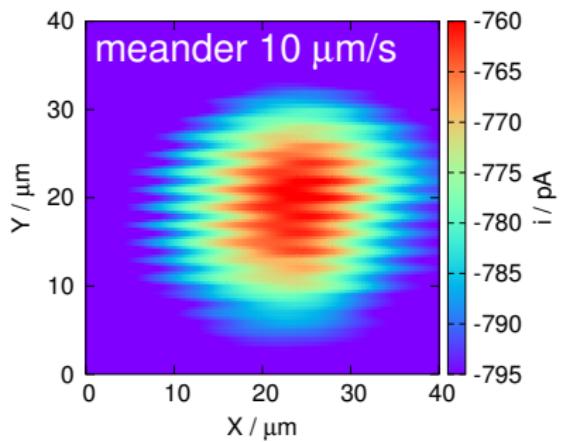

## Conclusion

1. SECM is a powerful tool. Many systems can be studied best with this technique.
2. However, it suffers from low scan speeds and high distortion.
3. Previously I have worked out a deconvolution technique to solve this problem for potentiometric SECM.
4. During my stay here I have worked out a similar method for amperometric SECM.
5. I have used the technique to restore images recorded by Monika and Phillip previously.
6. I have introduced single cell oxygen measurements here with the SECM, and successfully applied the deconvolution to those images as well.

# Thank you!

Dr. Monika Bozem  
Prof. Dr. Markus Hoth  
Phillip Knapp  
DAAD  
CIPMM

Thank you for your kind attention!